

# Varona Health Pillow **VIBE**

HK-7500

Varona Health Pillow, Vibe, induces sound sleep and promotes recovery from fatigue, helping the body recover its natural healing power.

 **HK Co.,Ltd.**

# Varona Health Pillow **VIBE** HK-7500



## ◎ 제품구성

베개, 아답타



## ◎ Product specifications and How to use

Items	Product specifications
Product name	Health pillow – Vibe
안전인증번호	HU10518-9003
정격전압	AC 100~240V, 50~60Hz
정격출력	DC5V
규격	390 x 130 x 70mm
중량	1.95Kg / 1EA
제조원	HK
제조국가	Made in Korea

## ◎ How to use

1. Plug the power cord of the main body of the product into a mains outlet.
2. Vibration begins at the same time when power is supplied.
3. The product continues to vibrate for 30 minutes and then comes to an automatic stop.
4. Disconnect the power plug and then connect it again if you want to use the pillow after the vibration is over.
5. Pull out the power plug when the use is over.




아답터 잭단자에 연결



# Varona Health Pillow VIBE

HK-7500

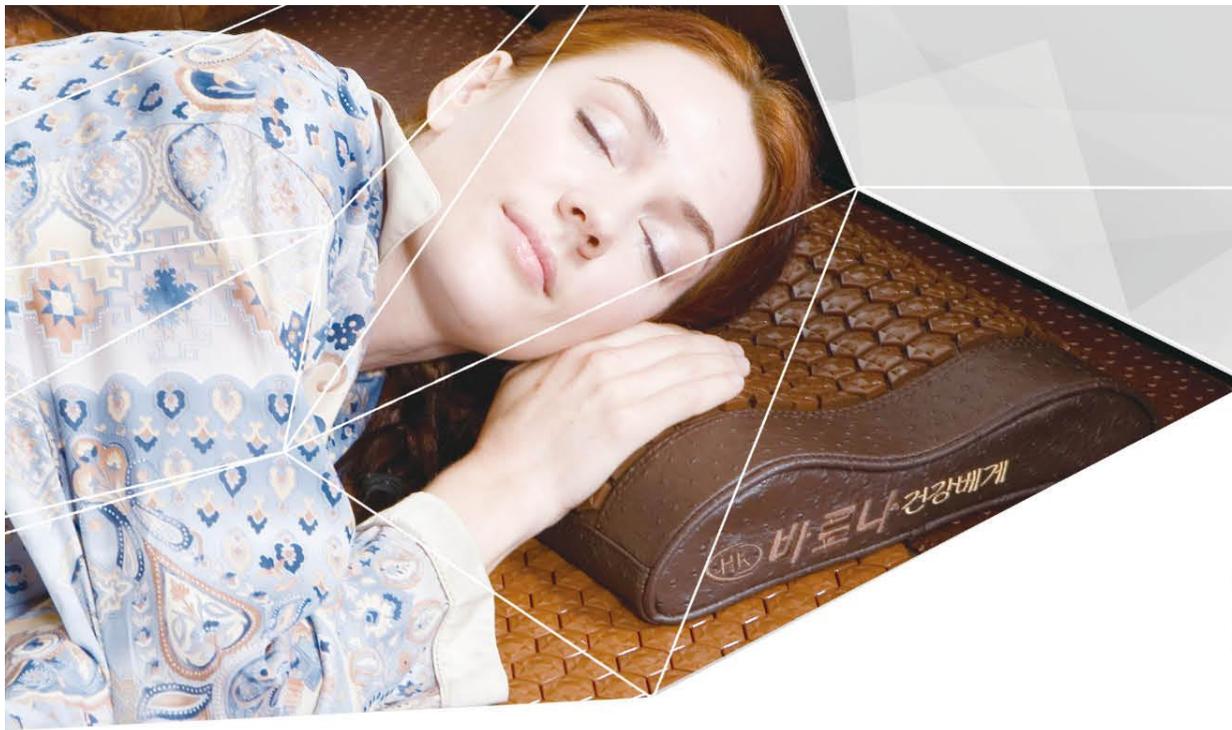

## ● Safety precautions

1. Check whether the power plug is disconnected if the product is to be cleaned.
2. Do not use wet objects near the product nor keep it in a moist place.
3. Do not prick or scratch the product with a sharp object like a knife or nail.
4. Do not disassemble the product at your discretion.
5. Disconnect the power plug when you leave your house or do not use the product for a long time.
6. Do not give strong impacts to the product or put heavy objects on it.

## ● Features of Varona Health Pillow Vibe

A C-shaped ergonomic design and excellent restoring capability

- Disperses pressure evenly across the cervical area to preserve its natural shape

Application of ceramic material produced of natural mineral ores using nano-technology and accumulation of heat through ceramic hexagons of a 3D pyramidal structure

- Emits a large amount of far infrared rays and anions continuously, promoting blood circulation

Soft vibration massaging

- Help the user have smooth blood circulation in the carotid artery and enjoy a sound sleep

Such features as ergonomic design to prevent the deformation of the cervical vertebrae, application of special ceramic material that continuously emits abundant far infrared rays and anions, and soft vibration massaging will ensure improved blood circulation and sound sleep for you.



## ● Ergonomic design

### Varona Health Pillow - Vibe


A C-shaped ergonomic design and the use of memory foam having excellent restoring power distribute the pressure inside the pillow so evenly as to prevent the neck bone from getting deformed while its soft vibration massaging helps the carotid artery flow smoothly and induces sound sleep.

### Regular pillow

Regular pillows that do not consider the shape of the cervical vertebrae or excessively cushioned pillows do not disperse the pressure inside evenly, thus resulting in the deformation of the neck bone or herniated cervical discs.

## ※ The importance of a pillow

- Regional symptoms resulting from pressure on cervical nerves

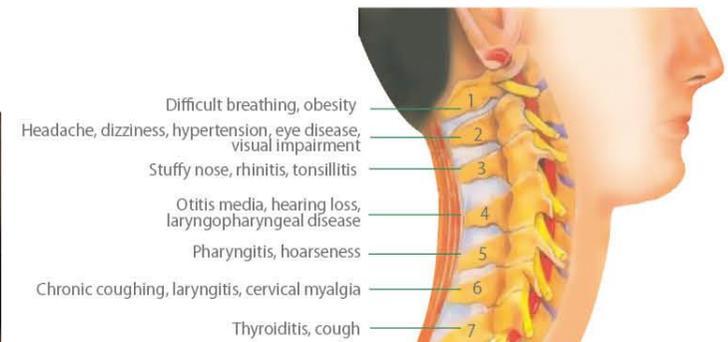

# Varona Health Pillow VIBE

HK-7500



## Effects of special ceramics

Energy-holding ceramics specially processed from natural ores by nano-technology

- Emits abundant far infrared rays and anions continuously.
- Promotes blood circulation and gets a sound sleep.

Heat-accumulating ceramics having a pyramid structure

- Maximizes far infrared effects through the embodiment of color changes according to temperatures



Biotite

Kiyoseki

Tourmaline

Germanium

Chunkwangsuk

kumkangyaktol

## Soft vibration massaging

- Soft vibration massaging felt during sleeping

## 원적외선

적외선 중 25 $\mu$ m 이상의 긴 파장으로 눈에 보이지 않고 열작용이 크며 침투력이 강하다. 특히, 유기화합물 분자에 대한 공진 및 공명 작용이 크다. 이러한 열작용은 각종 질병의 원인이 되는 세균을 없애고, 모세혈관을 확장시켜 혈액순환과 세포조직 생성에 도움을 준다. 또 세포를 1분에 2,000번씩 미세하게 흔들 어줌으로써 세포조직을 활성화하여 노화방지, 신진대사 촉진, 만성피로 등 각종 성인병 예방에 효과가 있다.

## 음이온

공기를 정화한다거나 인체의 자율신경계를 조절하는 효과가 있다고 하여 공기 중의 비타민이라고 지칭하기도 한다. 적당한 음이온 요법은 혈액을 중화하고 저항력을 증가시켜준다. 고온의 온열을 조사하여도 사용하면 시원함을 느끼는 것이 바로 음이온의 특징이다.

## Quality Warranty

Product Name : Varona Sofa Heating Pad

The above-mentioned product is warranted for its quality as described below:

Warranty period	One (1) year from the date of purchase
Date of purchase	YY    MM    DD
Customer	Address:
	Name:
	Tel. No.:
Store	Name of store:
	Name of director:
	Address
	Tel. No.:

1. This product undergoes our strict quality inspection procedure before it is shipped out.
2. This product is subject to free repair service for one year from the date of purchase according to the Consumer Protection Act if the defect is determined to result from normal use.
3. Contact the store or our A/S Center (1577-6778) or our Homepage (<http://www.hkmedical.co.kr>) if you have any problems with this product.
4. Services will be charged even during the warranty period, upon occurrence of any of the following events:
  - A. When the defect results from misuse or abuse or in connection with user's careless or negligent act.
  - B. When the defect results from willful alteration or repair by an authorized person.
  - C. When the product becomes defective or damaged as a result of fall.
  - D. When this warranty is not presented.
  - E. When this warranty does not contain the date of purchase and the seal and signature of the store where the product is purchased.
  - F. When the serial number marked on the inspection sticker has been damaged.

 **HK Co.,Ltd.**

### Manufacturer

Head office and factory :

669-2, Shinwol-ri, Chowol-eup, Gwangju-si, Gyeonggi-do, Korea

Tel. No. : 031-798-1577

Fax. No. : 031-798-1588

Contact Number for after-service : 1577- 6778

※ Issuing date: July 22, 2014 (Rev.0)